tribute®

social mobility

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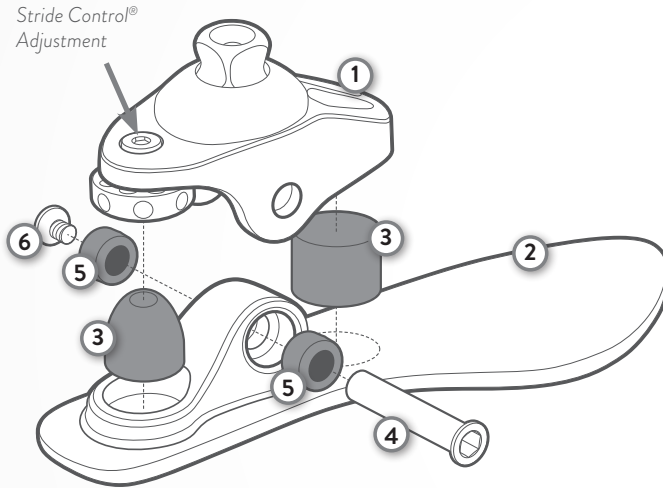
This diagram is to help familiarize you with the unique parts of the Tribute. These parts are referenced in the instructions and used when speaking with a technical service representative.

KEY COMPONENTS

1. Endo Ankle Bone
2. Foreheel Bone
3. Bumpers (2)
4. Axial Pin
5. Ankle Bushings (2)
6. Axial Pin Screw
Torque 4 N-m (36 in-lbs)

PACKAGE CONTENTS

- (1) Tribute Foot
- (1) Foot Shell
- (1) CPI Sock
- (1) Cosmetic Attachment Plate (regular toe only)
- (1) TruLube Lubricant

TOOLS REQUIRED

- (1) 4 mm Hex Key
- (1) 6 mm Hex Key

TOOLS RECOMMENDED

- Foot Horn
- Pin Guide

WARNING

- Do not expose this product to corrosive materials, salt water or pH extremes.
- Failure to follow these technical instructions or use of this product outside the scope of its Limited Warranty may result in injury to the patient or damage to the product.
- Any further disassembly or modification of components will void the warranty.

PROTECTIVE COVER ON DOME

Remove the protective cover on dome after alignment is completed and before patient leaves clinic.

ENDOSKELETAL MOUNTING

Use only high quality proximal endoskeletal components (30 mm).

ASSEMBLY AND DISASSEMBLY

1. Use the Foot Horn to don and doff the foot shell. Remove the CPI Sock and replace as needed.
2. Use a 4 mm hex key to reduce the preload on the bumpers by turning the Stride Control adjuster counter-clockwise. Record the number of turns for use in re-assembly.
3. Use a 6 mm and 4 mm hex keys to remove the axial pin screw.
4. Attach the pin guide to the axial pin. Push and rotate clockwise.
5. Unscrew the axial pin from the pin guide and remove it from the ankle bone.
6. Remove the ankle bone from the foreheel to access the front and rear bumpers and ankle bushings.
7. For re-assembly, lubricate the axial pin, inside and outside of the ankle bushings and ankle bushing pocket, then reverse steps 1-4. Torque axial pin screw to 4 N·m (36 in·lbs).

DO NOT LUBRICATE FRONT AND REAR BUMPERS.

STATIC ALIGNMENT

For optimal function have the patient's weight balanced between the forefoot and heel. The load line divides the foot at 1/3 heel lever and 2/3 toe lever.

STRIDE CONTROL ADJUSTMENT

The adjustable Stride Control gives you the ability to customize an individual's gait timing, fine-tuning the foot's response by either turning it clockwise/ counterclockwise. This adjustment affects both plantar and dorsiflexion resistance without the need to change soft components.

The Tribute ships with a medium Stride Control setting. If new bumpers are fitted, set the preload to 2-turns (medium preload) from the point of initial contact with the ankle bone.

If an additional 1-1/2 turns are required to attain the desired resistance, then a change in bumpers is recommended.

DYNAMIC ADJUSTMENTS

DESIRED RESULT	ALIGNMENT CHANGE	COMPONENT CHANGE
Firmer Toe Response	Plantarflex the Tribute or move load line posterior	Increase Stride Control Front Bumper up one step firmer
Softer Toe Response	Dorsiflex the Tribute or move load line anterior	Decrease Stride Control Front Bumper down one step softer
Firmer Heel Response	Dorsiflex the Tribute or move load line anterior	Increase Stride Control Rear Bumper up one step firmer
Softer Heel Response	Plantarflex the Tribute or move load line posterior	Decrease Stride Control Rear Bumper down one step softer

GAIT MATCHING® GUIDELINES

Contact College Park Technical Service if you have questions about modifying the gait match.

WARRANTY INSPECTION / MAINTENANCE INFORMATION

College Park recommends that you schedule your patients for check-ups on the Tribute foot six months after initial fitting and then annually.

High patient weight and/or impact level may require more frequent inspections. Soft component wear depends on the patient weight, impact level and environment. We recommend you inspect the following applicable parts for excessive wear and fatigue at each warranty inspection and replace as needed.

- Soft components (disassemble, inspect & re-lubricate)
- Composites and Adapters
- CPI Sock
- Foot Shell

TECHNICAL ASSISTANCE / EMERGENCY SERVICE 24-7-365


College Park's regular office hours are Monday through Friday, 8:30 am – 5:30 pm (EST). After hours, an emergency Technical Service number is available to contact a College Park representative.

LIABILITY

The manufacturer is not liable for damage caused by component combinations that were not authorized by the manufacturer.

CAUTION

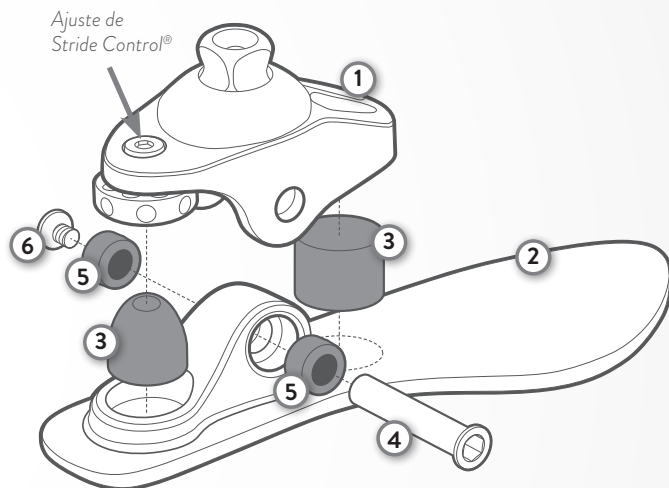College Park products and components are designed and tested according to the applicable official standards or an in-house defined standard when no official standard applies. Compatibility and compliance with these standard is achieved only when College Park products are used with other recommended College Park components. This product has been designed and tested based on single patient usage. This device should NOT be used by multiple patients. If any problems occur with the use of this product, immediately contact your medical professional.

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©2018 College Park Industries, Inc. All rights reserved. College Park, Stride Control, Gait Matching
and Tribute are registered trademarks of College Park Industries.