

TPE Liner Family

Experience the difference!



Nothing comes closer to your patient than their liner



The liner defines the feeling a person has for his or her prosthesis. The liner needs to do more than feel good. It needs to protect and cushion their limb as well from the first use and as long as they wear the liner.

We've drawn on deep experience from technicians around the world to improve the TPE-liner to better meet the needs of your patients. The result is a new gel formulation that

- is highly elastic
- makes it easier for patients to flex their knee
- contains white oil that help relieve dry skin
- a TPE that can be delivered to patients with unusual limb shapes
- is available in locking and cushion versions



6Y90/6Y92 TPE liner

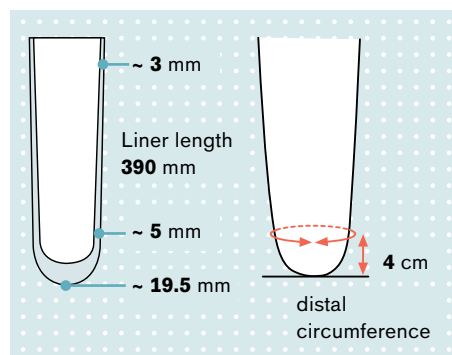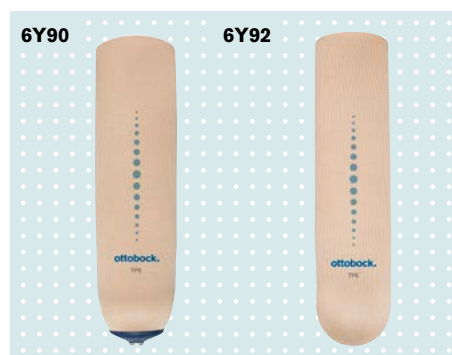

Article no.	6Y90	6Y92
Amputation Level	transtibial (TT)	transtibial (TT)
Material	Copolymer (TPE)	Copolymer (TPE)
Distal connection	With	Without
Textile cover	With	With
Socket design	Total surface weight-bearing socket	Total surface weight-bearing socket
Matrix	10 cm length	Without
Sizes (see distal circumference)	200 mm, 250 mm, 280 mm, 320 mm, 260 mm	

i Measure the residual limb circumference 4 cm above the distal end of the residual limb and select the liner size which corresponds to this measurement, or the next smallest size.

Order no. (ordering example):


Reference number	=	Size
6Y92	=	280

Characteristics and advantages of 6Y90/92 TPE Liners:

- Transtibial liner
- 3mm tapered liner
- The 6Y90 TPE liner contains an integrated matrix in the distal area. This helps prevent the so-called “milking effect.”

6Y93 Balance liner

6Y93

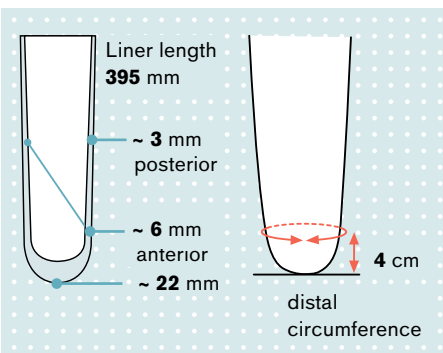

Article no.	6Y93=L6			6Y93=C6	
Balance	Thicker front side with 6 mm wall thickness for enhanced user comfort, thinner posterior side with 3 mm wall thickness to facilitate flexion				
Amputation Level	Transtibial (TT)			Transtibial (TT)	
Material	Copolymer (TPE)			Copolymer (TPE)	
Distal connection	With			Without	
Textile cover	With			With	
Socket design	Total surface weight-bearing socket				
Matrix	Without			Without	
Sizes	4 cm from distal			30 cm from distal (reference)	
		Min	Max	Min	Max
	S	150	220	200	270
	M	180	260	230	360
	MP	200	280	330	500
	L	250	330	360	550
	LP	280	400	400	600
	XL	330	450	430	650

i Measure the residual limb circumference 4 cm above the distal end of the residual limb and select the liner size which corresponds to this measurement, or the next smallest size.

Order no. (ordering example):

Reference number	=	Size
6Y93	=	C6-M

Characteristics of the 6Y93 Balance TPE liners:

- Transtibial liner
- 6 mm tapered liner with thicker anterior for protection and thinner posterior for flexibility

Ottobock
U.S. P 800 328 4058 · F 800 962 2549
Canada P 800 665 3327 · F 800 463 3659
www.ottobockus.com · www.ottobock.ca