



PDAC verified L1843/L1851*

GenuTrain® OA

Orthosis for compression and unloading of the knee joint.

GenuTrain OA with its innovative unloading system, relieves the lateral or medial compartment of the knee based on the 3-point principle. The unloading relief is easily adjustable with the Boa Fit-System. A low-profile and light weight design, combined with 3D AirKnit compression support, allows for unparalleled relief, joint control and all day wearing comfort.

Indications

- Medial or lateral osteoarthritis of the knee
- Instability of the knee joint (mediolateral)
- Unicompartamental cartilage/bone lesion
- Condition after cartilage surgery and cartilage growth for necessary unilateral post-operative relief
- Following unilateral meniscal lesions
- Ahlbäck's disease (osteonecrosis)
- Osteochondrosis dissecans
- Bone bruising (bone edema of the medial condylus of the femur/tibial plateau)

*Suggested codes are based on approved coding from PDAC as of 01/20. It is the responsibility of the provider to determine the appropriate billing code, as well as, whether the use of a product complies with medical necessity and other documentation requirements of the payer.

The new GenuTrain® OA

MORE ACTIVITY, LESS PAIN.



SUCCESSFUL THERAPY DRIVEN BY HIGH PATIENT COMPLIANCE

TOOL-FREE:

- Extension 0°, 15°, 30°, 45°, 60°
- Flexion 30°, 45°, 60°, 75°, 90°




INNOVATIVE UNLOADING SYSTEM

- Adjustable power-8 strap system and stabilizing hinge provide relief (3 point principle)



TOOLFEE ADJUSTMENT

- Hinge automatically adjusts to leg shape (varus/valgus) and allows for flexion and extension control

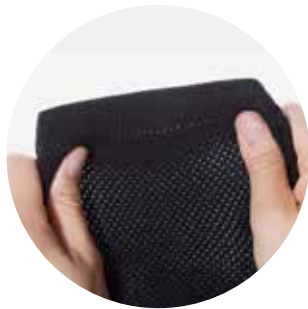

COMPRESSION AIR KNIT

- Compression knit for improved proprioception
- Extra soft popliteal area
- Anatomically contoured for perfect fit
- Latex-free



EXCELLENT WEARING COMFORT

- Lightweight, low profile design
- Extra breathable 3D Airknit
- Strap free popliteal area for pressure free bending zone
- Easy and secure donning and doffing with ergonomic finger loops




BOA® FIT-SYSTEM

- Easy adjustable unloading relief during activity



Limited inventory required due to design allowing for medial and lateral OA relief

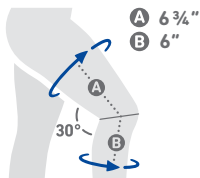

Size	1	2	3	4	5	Order number
Circumfer. 6 3/4" above knee	15 3/4 - 18 1/2	18 1/2 - 20 1/2	20 1/2 - 22 3/4	22 1/2 - 25 1/2	23 1/2 - 28 1/4	12141650070+0°
Circumfer. 17 cm above knee	40 - 47	47 - 52	52 - 58	57 - 65	60 - 72	
Circumfer. 6" below knee	11 - 13 1/2	13 1/2 - 15 1/4	15 1/4 - 17	17 - 18	18 - 19 3/4	
Circumfer. 15 cm below knee	28 - 34	34 - 39	39 - 43	43 - 46	46 - 50	

+6 = right, +7 = left, ° = size

BAUERFEIND USA, INC. 1230 Peachtree St. NE • Suite 1000 • Atlanta, GA 30309
P 1 800 423 3405 • P (404) 201 7800 • F (404) 201 7839
E-mail info@bauerfeindusa.com