

# Guardian

GLAK(HD) / GLBK(HD)



Transfemoral (AK) rings begin 5 1/2 inches from the center of the distal end.

Transtibial (BK) rings begin 3 1/2 inches from the center of the distal end.

## *Gel Liner, Reinforced with Gel Bands – No Distal Attachment*

The ALPS™ Guardian Suction Liner features raised bands that grip the socket wall to form a secure interface between the socket and the liner. Low modulus bands stretch against the socket wall, while the smooth inner wall conforms easily to the residual limb, to ensure there is no restriction of blood flow or stiffening to inhibit donning.

- Available with GripGel and HD Gel
- Full length inversion for ease of donning
- No reinforced seams, single piece construction that still ensures durability
- Available in 3mm and 6mm uniform



Cushion



GripGel



High Density (HD) Gel



Black Fabric,  
Horizontal and  
Vertical Limited Stretch

Use the Skin Reliever (ENCP) with the Guardian Liner to compensate for volume fluctuations with continued use of same socket.




US Patent No. 8,852,291



Tel: 727.528.8566 | Tel: 800.574.5426  
www.easyliner.com | info@easyliner.com

## Cross Section Diagram

3mm Uniform, 6 Uniform

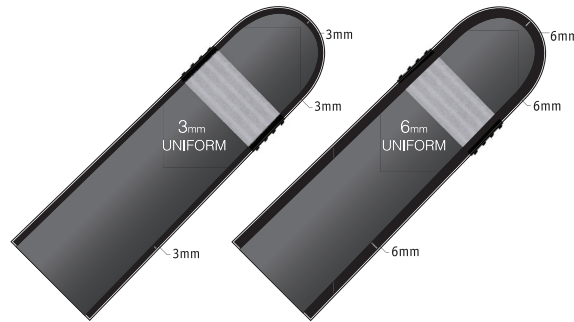

## How to Size


Measurements taken 6cm from the distal end

### Transfemoral Guardian Suction Cushion Gel Liner

Measured Size (cm) Description	GLAK 3mm Uniform	GLAK 6mm Uniform
26-27 cm	GLAK(HD) <b>26-3</b>	GLAK(HD) <b>26-6</b>
28-31 cm	GLAK(HD) <b>28-3</b>	GLAK(HD) <b>28-6</b>
32-37 cm	GLAK(HD) <b>32-3</b>	GLAK(HD) <b>32-6</b>
38-43 cm	GLAK(HD) <b>38-3</b>	GLAK(HD) <b>38-6</b>

### Transtibial Guardian Suction Cushion Gel Liner

Measured Size (cm) Description	GLBK 3mm Uniform	GLBK 6mm Uniform
16-19 cm	GLBK(HD) <b>16-3</b>	GLBK(HD) <b>16-6</b>
20-23 cm	GLBK(HD) <b>20-3</b>	GLBK(HD) <b>20-6</b>
24-25 cm	GLBK(HD) <b>24-3</b>	GLBK(HD) <b>24-6</b>
26-27 cm	GLBK(HD) <b>26-3</b>	GLBK(HD) <b>26-6</b>
28-31 cm	GLBK(HD) <b>28-3</b>	GLBK(HD) <b>28-6</b>

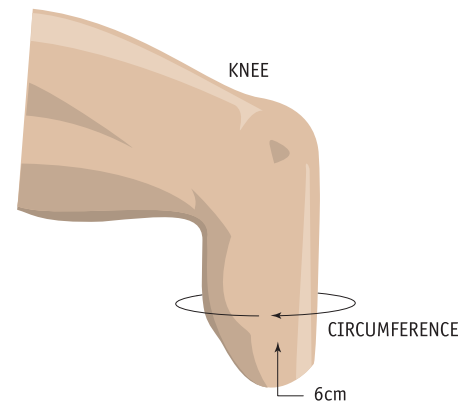

## Skin Reliever (ENCP) Sizing Chart

Measurements taken 6cm from the distal end



Use the Skin Reliever (ENCP) with the Guardian Liner to compensate for volume fluctuations with continued use of same socket.

Measured Size (cm) Description	ENCP Tapered 3mm-1.5mm
10-19 cm	ENCP- <b>16</b>
20-27 cm	ENCP- <b>20</b>
28-37 cm	ENCP- <b>28</b>