Maverick junior. A smile with every step.

The F24 Maverick junior was designed for active children and adolescents who navigate varied indoor and outdoor environments and place a high value on flexibility and durability.

Designed for optimal outcomes

- The fiberglass design delivers a high degree of durability and flexibility with the expected energy return
- The universal rounded toe element offers users increased robustness
- Low build height
- Pyramid or lamination anchor connection options
- The integrated heel wedge with optional insert makes it possible to customize the shock absorption at heel strike
- The sole of the footshell features a smiley face as well as a texture for walking barefoot
- Easy cleaning
- Waterproof Protected against fresh, salt, and chlorinated water

More than a foot. A foundation.





Technical data

For	Children, adolescents and adults with small feet
Max. body weight	99 lbs (45 kg)
Sizes	16 – 21 cm
Weight without footshell*	7.4 oz (210 g)
Footshell shape	Slim shape for a heel height of 5 mm
Footshell color	beige (4), brown (15)
Weight w/footshell*	12.4 oz (352 g)
Build height*	1 7/8" (47 mm)

^{*} Reference size is 21 cm.



F24 *Maverick junior*. Order form.


Customer no.		Customer no.		
Company		Company		
Street		Street		
Postal/ Zip code	City	Postal/ Zip code	City	
O&P professional		Customer ref.		

Selection of the spring stiffness relative to body weight and foot size

Foot size Body weight	16 cm	17 cm	18 cm	19 cm	20 cm	21 cm
20 - 34 lbs (9-15 kg)	1	1	1	1	-	-
35 - 53 lbs (16-24 kg)	2	2	2	2	2	-
54 - 75 lbs (25-34 kg)	-	3	3	3	3	3
76 - 99 lbs (35-45 kg)	_	-	4*	4*	4*	4*

^{*} Stiffness 4 is not available for the Maverick junior with lamination adapter.

Example for ordering F24=L-16-2-P/4N

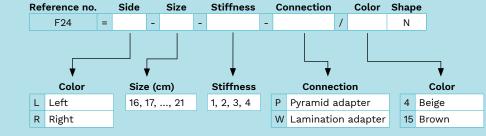



Scope of delivery

The scope of delivery includes the F24 foot module, the 2C24 footshell, the 2F14 heel wedge insert, the 2C101 footshell removal tool, a black Spectra sock as well as a quick guide for users and their families on how to remove the footshell for easy cleaning.

For the F24 foot module with lamination adapter, the 2R89=1 lamination anchor and the 601T24=* cap screw are also included in the scope of delivery.

Recommended combinations **

TT	TF
	3R67
Children modular system 4R11=* Quickchange 2R41=1 4R66 4R110 4R60	4R11=* Quickchange 2R41=1 4R66 4R110 4R60
Adult modular system (waterproof) 4R11=* Quickchange 2R37 with 4X907=1 Pediatric to Adult Adapter 2R38 with 4X907=1 Pediatric to Adult Adapter	

^{**} All components are sold separately and are available Ottobock products that are compatible with the F24 Maverick junior, which helps ensure optimal performance. Practitioners need to select components based upon individual patient criteria.

The Maverick junior prosthetic foot is also approved for combination with Genium X4, C-Leg and Kenevo. The selection of the additional prosthetic parts must be matched to the knee joints with a kids to adult adapter by the technician.