



**Made to
move.**

Dynamic Walk

Fillauer®

Reimagining Possibilities.

The Dynamic Walk AFO is designed to allow people with dorsiflexor paresis or paralysis to confidently regain their stride. Whether their goal is to walk, hike, jog, or play, the Dynamic Walk's flexible dorsi-assist joints allow them to move freely.

Easy Alterations

Both the cuff and the foot plate are made of a lightweight thermocomposite that is easily trimmable and remoldable

Better Movement

Flexible PEEK rods are available in a variety of configurations for the perfect amount of dorsi-assist

Dynamic Foot Plate

The thermocomposite footplate is flexible for multi-planar movement and easily fits in most shoes





Standard



Wide Ankle



Double-Double

Dynamic Walk, Standard

Size	Length	Width	Height	Calf Circumference	Weight	Left	Right
Small	9.1 in. (230 mm)	2.6 in. (67 mm)*	12.6 in. (320 mm)	9.4 – 13.0 in. (240 – 330 mm)	3.7 oz. (105 g)	700111121	700111122
Medium	9.8 in. (250 mm)	2.8 in. (72 mm)*	14.0 in. (355 mm)	13.4 – 16.5 in. (340 – 420 mm)	4.3 oz. (122 g)	700111141	700111142
Large	10.6 in. (270 mm)	3.1 in. (78 mm)*	14.6 in. (370 mm)	16.9 – 20.5 in. (430 – 520 mm)	5.1 oz. (144 g)	700111161	700111162

Dynamic Walk, Wide Ankle

Size	Length	Width	Height	Calf Circumference	Weight	Left	Right
Small	9.1 in. (230 mm)	2.8 in. (72 mm)*	13.2 in. (335 mm)	9.4 – 13.0 in. (240 – 330 mm)	3.7 oz. (106 g)	700111321	700111322
Medium	9.8 in. (250 mm)	3.1 in. (80 mm)*	14.0 in. (355 mm)	13.4 – 16.5 in. (340 – 420 mm)	4.3 oz. (123 g)	700111341	700111342
Large	10.6 in. (270 mm)	3.4 in. (87 mm)*	15.0 in. (380 mm)	16.9 – 20.5 in. (430 – 520 mm)	5.0 oz. (142 g)	700111361	700111362

Dynamic Walk, Double-Double

Size	Length	Width	Height	Calf Circumference	Weight	Left	Right
Small	9.1 in. (230 mm)	2.7 in. (68 mm)*	12.6 in. (320 mm)	9.4 – 13.0 in. (240 – 330 mm)	5.5 oz. (157 g)	700117121	700117122
Medium	9.8 in. (250 mm)	2.9 in. (73 mm)*	14.6 in. (370 mm)	13.4 – 16.5 in. (340 – 420 mm)	6.6 oz. (188 g)	700117141	700117142
Large	10.6 in. (270 mm)	3.1 in. (80 mm)*	15.4 in. (400 mm)	16.9 – 20.5 in. (430 – 520 mm)	7.5 oz. (213 g)	700117161	700117162

*Width measured approximately one centimeter above the foot plate

Padding Kit

When ordering, please note that the padding must be ordered separately.

Size	Product Number
Small	700100220
Medium	700100240
Large	700100260





Single Side Medial



Single Side Lateral



Shoe Mount

Dynamic Walk, Single Side Medial

Size	Length	Height	Calf Circumference	Weight	Left	Right
Small	9.1 in. (230 mm)	13.4 in. (340 mm)	9.4 – 13.0 in. (240 – 330 mm)	3.6 oz. (103 g)	700113121	700113122
Medium	9.8 in. (250 mm)	14.6 in. (370 mm)	13.4 – 16.5 in. (340 – 420 mm)	4.2 oz. (120 g)	700113141	700113142
Large	10.6 in. (270 mm)	15.2 in. (385 mm)	16.9 – 20.5 in. (430 – 520 mm)	4.9 oz. (140 g)	700113161	700113162

Dynamic Walk, Single Side Lateral

Size	Length	Height	Calf Circumference	Weight	Left	Right
Small	9.1 in. (230 mm)	13.4 in. (340 mm)	9.4 – 13.0 in. (240 – 330 mm)	3.6 oz. (103 g)	700114121	700114122
Medium	9.8 in. (250 mm)	14.6 in. (370 mm)	13.4 – 16.5 in. (340 – 420 mm)	4.2 oz. (120 g)	700114141	700114142
Large	10.6 in. (270 mm)	15.2 in. (385 mm)	16.9 – 20.5 in. (430 – 520 mm)	4.9 oz. (140 g)	700114161	700114162

Dynamic Walk, Shoe Mount

Size	Length	Height	Calf Circumference	Weight	Left	Right
Small	6.5 in. (165 mm)	13.4 in. (340 mm)	9.4 – 13.0 in. (240 – 330 mm)	3.5 oz. (100 g)	700110721	700110722
Medium	6.5 in. (165 mm)	14.6 in. (370 mm)	13.4 – 16.5 in. (340 – 420 mm)	3.9 oz. (110 g)	700110741	700111142
Large	6.5 in. (165 mm)	15.2 in. (385 mm)	16.9 – 20.5 in. (430 – 520 mm)	4.2 oz. (120 g)	700110761	700110762

Padding Kits

When ordering, please note that the padding must be ordered separately.

Size	Medial Left	Medial Right	Lateral Left	Lateral Right
Small	700100322	700100321	700100321	700100322
Medium	700100342	700100341	700100341	700100342
Large	700100362	700100361	700100361	700100362





Dynamic Walk Custom

The **Custom Dynamic Walk** is recommended for patients who need to accommodate the unique shape of their foot and leg and/or specific functional needs. Whether requiring greater toe lift or unique sizing, including pediatric sizing, Fillauer's Custom Dynamic Walk will be created to suit your patients' clinical needs.

Custom fabrication forms are available on the following pages. For more information, contact a Fillauer Inside Sales Specialist at **423.624.0946** or **centralfab@fillauer.com**.

Dynamic Walk: Adult

Custom Fabrication Order Form

Fillauer®

2710 Amnicola Highway
Chattanooga, TN 37406
423.624.0946
Fax 423.698.6076
centralfab@fillauer.com

Company

Practitioner

Account No.

Phone

Address Line 1

Email

Address Line 2

Patient ID

City, State, Zip

Purchase Order

Diagnosis

700100100 Custom Dynamic Walk

☐ Standard (2 PEEK Rods)	☐ Double Double (4 PEEK Rods)
☐ Single Side Lateral ☐ Single Side Medial (2 PEEK Rods)	☐ Attach to Shoe ☐ Standard ☐ Medial Only ☐ Lateral Only

Please specify:

☐ Right Foot ☐ Left Foot ☐ Bilateral

Choose Options

- ☐ Carbon standard layup: same as custom fit version
- ☐ Reinforced foot plate for increased durability
- ☐ Reinforced calf band for increased durability
- ☐ Custom foot orthotic (+ \$40)
- ☐ Additional PEEK rods (+ \$50)

☐ Would you like us to contact you prior to fabricating the brace?

Notes _____


Mold Type

Negative cast, insole tracing, and measurement form preferred.

- ☐ Positive cast
- ☐ Negative cast
- ☐ Modified positive cast
- ☐ 3D Scan (.stl, .xyz, .obj, .op3 file types accepted)

Please note that if non-standard trim lines are not indicated, then standard trim lines will be applied.

Please fill in all of the measurements below

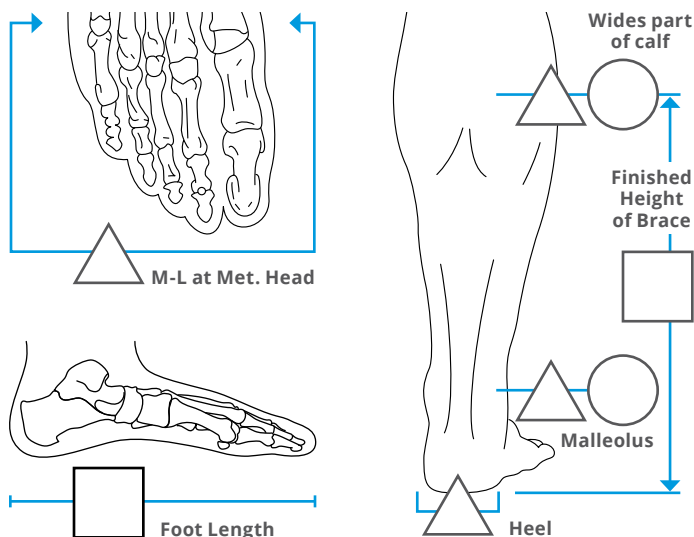

Casting Instructions

- **Include a tracing of the Insole of the shoe to be used.**
- **Please mark the desired height of brace on the finished cast.**
- Patient should be weight bearing, standing straight with upright posture when casting. The calf needs to be at a 90° angle to the horizontal plane with the foot at the required heel height. The patient needs to be casted up to the tibial tuberosity.
- Please mark the pressure points or any areas we should avoid covering on the positive or negative cast.

Dynamic Walk: Pediatric

Custom Fabrication Order Form

Fillauer®

2710 Amnicola Highway
Chattanooga, TN 37406
423.624.0946
Fax 423.698.6076
centralfab@fillauer.com

Company _____

Practitioner _____

Account No. _____

Phone _____

Address Line 1 _____

Email _____

Address Line 2 _____

Patient ID _____

City, State, Zip _____

Purchase Order _____

Diagnosis _____

700100100 Custom Dynamic Walk

☐ Standard (2 PEEK Rods)	☐ Double Double (4 PEEK Rods)
☐ Single Side Lateral ☐ Single Side Medial (2 PEEK Rods)	

Please specify:

☐ Right Foot ☐ Left Foot ☐ Bilateral

Choose Options

- ☐ Carbon standard layup: same as custom fit version
- ☐ Reinforced foot plate for increased durability
- ☐ Reinforced calf band for increased durability
- ☐ Custom foot orthotic (+ \$40)
- ☐ Additional PEEK rods (+ \$50)
- ☐ SMO (+ \$75)

☐ Would you like us to contact you prior to fabricating the brace?

Notes _____


Mold Type

Negative cast, insole tracing, and measurement form preferred.

- ☐ Positive cast
- ☐ Negative cast
- ☐ Modified positive cast
- ☐ 3D Scan (.stl, .xyz, .obj, .op3 file types accepted)

Please note that if non-standard trim lines are not indicated, then standard trim lines will be applied.

Please fill in all of the measurements below

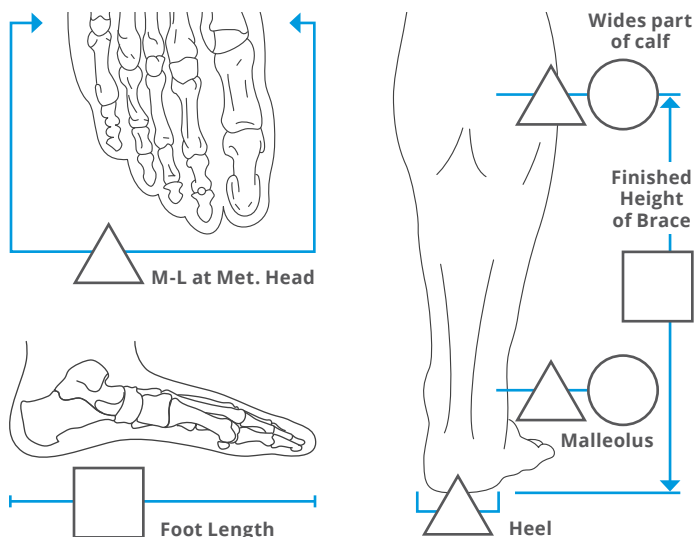

Casting Instructions

- **Include a tracing of the Insole of the shoe to be used.**
- **Please mark the desired height of brace on the finished cast.**
- Patient should be weight bearing, standing straight with upright posture when casting. The calf needs to be at a 90° angle to the horizontal plane with the foot at the required heel height. The patient needs to be casted up to the tibial tuberosity.
- Please mark the pressure points or any areas we should avoid covering on the positive or negative cast.



www.fillauer.com

Fillauer LLC

2710 Amnicola Highway
Chattanooga, TN 37406
423.624.0946
info@fillauer.com

© 2021 Fillauer LLC
AD364/01-28-16/02-01-21