



THUASNE

SafeLimb

Transtibial Residual Limb Protector



Features

Post-surgical
support & protection

Five sizes; universal
for left or right leg

Bars can be
hand-contoured &
lengthened

Anterior opening
& straps for easy
application

Cushioned rigid plate
under residual limb

SafeLimb


Transtibial Residual Limb Protector

RECOMMENDED L5450


During post-surgical care of transtibial amputation patients, the SafeLimb provides supportive protection of the residual limb if the patient falls.

Prosthetists will appreciate the SafeLimb's unique design, materials and adaptability – which enhances the fit, suspension and patient compliance. A spherical compressive foam cushion mounted on a rigid plate – attached to the struts at the distal end of the brace – prevents trauma to the surgical wound.

The length of the SafeLimb's medial and lateral bars is adjustable, and the bars can also be hand-contoured to align with the degree of extension. The material interface is lightweight and breathable, with anterior-opening and anterior straps for ease of application. A proprietary "SupraPatellar Strap" (above the knee) effectively prevents migration.

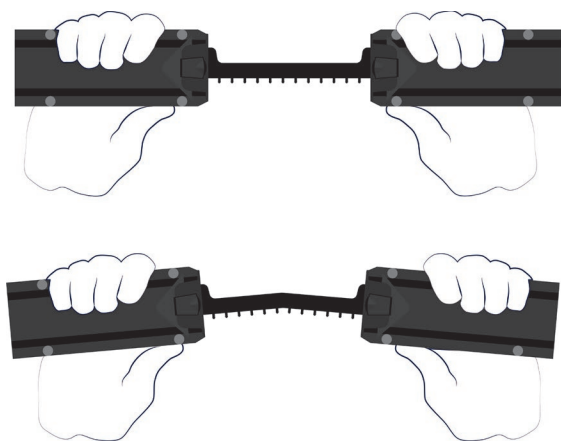

The SafeLimb's metal bars can be contoured without tools.

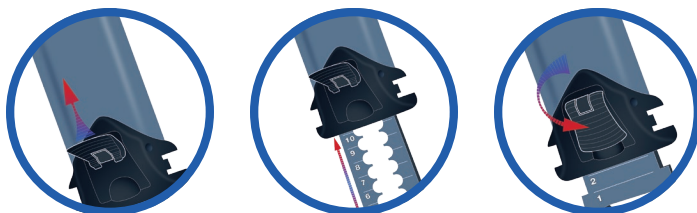

Townsend's proprietary "SnapLock" allows for quick, tool-free length adjustment.



The residual limb is protected by a cushioned plate at the distal end of the brace

Size	Thigh Circumference 6" Above Mid Patella
SM	14 – 18 ³ / ₄ "
MD	18 ³ / ₄ – 20 ³ / ₄ "
LG	20 ³ / ₄ – 24 ³ / ₄ "
XL	24 ³ / ₄ – 29 ¹ / ₂ "
XXL	29 – 34"



 **TOWNSEND**
THUASNE USA

Visit us at ThuasneUSA.com